

THE WIRELESS VIDEO OTOSCOPE



Delfino

Delfino

Delfino is the evolution of the traditional otoscope: the small and ergonomic hand piece allows the operator to see, capture and store top-quality images of the ear canal. All this with the maximum patient and operator comfort, thanks to the wireless feature.




Delfino is composed of the receiver base and the hand piece.

DESIGNED, ENGINEERED AND
MANUFACTURED IN ITALY



THE HAND PIECE

- Endoscopic gradient-index optics
- LED illumination
- CMOS camera
- Li-Ion battery, allowing two hours of uninterrupted use
- Wireless video transmission system
- Compatible with standard speculum.

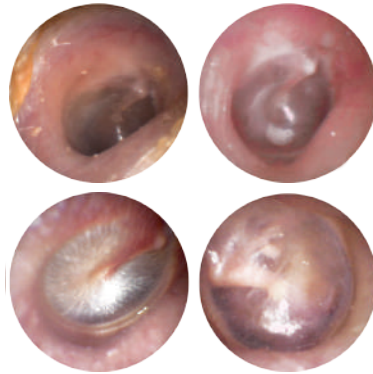

THE RECEIVER BASE

- Recharges the hand piece battery
- Embeds a pre-heating system for preventing the optics fogging
- USB digital output for the direct connection to the computer
- Analog video output for the connection to any standard display.

NOAH & MAESTRO INTEGRATION

Delfino comes with Maestro software, through which you can capture images and store them in the patient folder, together with other ENT and audiological examinations. If you are a Noah user, Maestro can be also installed as a Noah module.



TECHNICAL DATA

Illumination.....	LED array
Optics.....	Endoscopic gradient-index optics
Camera.....	CMOS sensor, res. 720 x 576
Min. working dist.....	10 mm
Battery.....	Li-Ion rechargeable
Video transmission.....	7 selectable channels
Pre-heating.....	in the receiver base
Video outputs.....	digital USB 2.0 and analog composite
Weight.....	base: 500 g - hand piece: 220 g
Dim. (WxHxD).....	base: 105 mm x 75 mm x 250 mm hand piece: 220 mm x 90 mm x 140 mm
Classification.....	Class I (MDD 93/42, Annex IX)

COMPUTER MINIMUM REQUIREMENTS

Processor.....	Intel i3 or similar
RAM.....	2 GB
Hard disk.....	1 GB free space
Display resolution.....	1024 x 768
Operating system.....	Microsoft Windows 7, 8 or 10



Delfino is developed by:

INVENTIS s.r.l.
CORSO STATI UNITI, 1/3
35127 PADOVA - ITALY
PHONE: +39.049.8962 844
FAX: +39.049.8966 343
info@inventis.it
www.inventis.it

Follow us on Twitter
<http://twitter.com/inventissrl>



The Inventis Quality System complies with ISO 9001 and ISO 13485 standards.

Inventis® is a registered trademark of Inventis s.r.l.

Australian Distributor:
SONIC Free Call:
Everyday Sounds Better 1800 639 263